

## Pro Series Class II B2 Biological Safety Cabinet



**LCD Display**  
7 inch touch screen display is easy to monitor all the safety parameters at a glance and ergonomically sized control panel improves user interface.



**Worktop Holder**



**Worktop Holder Usage Diagram**




**Work Zone**  
Work zone, made of 304 stainless steel, is surrounded by negative pressure.



**Remote Control**  
All functions can be realized with it, making the operation much easier and more convenient.



**UV Lamp**  
Emission of 253.7 nanometers for most efficient decontamination.

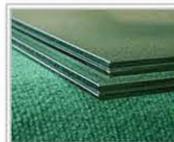

**Front Window**  
Two-layer laminated toughened glass ≥ 6mm, Anti UV



**Footmaster Caster**  
Universal caster with brake and leveling feet



**Foot Switch**  
Adjust front window height by foot during experiment, to avoid airflow turbulence caused by arm movement.



**Waterproof Socket**  
2 waterproof sockets are located in the side wall, for optimum convenience of using small devices inside the cabinet.



## ADVANTAGE

1. Larger angle of view,easy to operate and more humanized.
2. Motorized front window.
3. ULPA filter life and UV life indicator.
4. Time reserve function.
5. With memory function in case of power-failure.
6. Self-locking button switch to solve the trouble of losing the key.
7. 7 inch Large touch screen with Intelligent control system: operators can check detailed status of the cabinet, such as inflow and downflow velocity, work area temperature and humidity, filter pressure, UV working time and filter working time, filter life indicator etc.
8. Audio and visual alarm: filter replacement, window over height, abnormal air flow velocity, etc.High filter pressure alarm.
9. Remote control: all functions can be realized with it, making the operation much easier and more convenient.
10. Interlock function: UV lamp and front window; UV lamp and blower, LED lamp; blower and front window
11. Foot switch: adjust front window height by foot during experiment, to avoid airflow turbulence caused by arm movement.

## TECHNICAL PARAMETERS

Model	BSC-1100IIB2-Pro	BSC-1500IIB2-Pro	BSC-1800IIB2-Pro	
External Size (W*D*H)	1100*750*2255mm	1500*770*2290mm	1875*795*2280mm	
Internal Size (W*D*H)	940*600*655mm	1350*600*660mm	1700*600*660mm	
Tested Opening	Safety Height 200mm(8")			
Max Opening	420 mm (17")	500mm (20")	480mm (20")	
Inflow Velocity	0.53±0.025m/s			
Downflow Velocity	0.33±0.025m/s			
ULPA Filter	Two, 99.9995% efficiency at 0.12µm, filter life indicator.			
Front Window	Motorized, two-layer laminated toughened glass ≥ 6mm, anti UV			
Noise	≤ 67dB			
UV Lamp	30W*1	40W*1	40W*1	
	UV timer, UV life indicator, emission of 253.7 nanometers for most efficient decontamination.			
Illuminating Lamp	LED lamp			
	12W*2	16W*2	16W*2	
Illumination	≥1000Lux			
Consumption	800W	1300W	1400W	
Waterproof Sockets	Two, total load of two sockets: 500W			
Display	7 inch touch screen display: exhaust filter and downflow filter pressure, filter and UV lamp working time, inflow and downflow velocity, filter life, humidity and temperature, system working time etc.			
Control System	Microprocessor			
Airflow System	0% air re-circulation, 100% air exhaust			
Alarm	Abnormal airflow velocity; Filter replacement; Front window at unsafe height; High filter pressure alarm.			
Exhaust Duct	4 meters PVC duct, Diameter:300mm			
Display	Work Zone: 304 stainless steel Main Body: Cold-rolled steel with anti-bacteria powder coating			
Work Surface Height	750mm(Size can be customized according to requirements)			
Caster	Footmaster caster			
Power Supply	AC 220V±10%, 50/60Hz; 110V±10%, 60Hz			
Standard Accessory	LED lamp*2, UV lamp*2, base stand, remote control, foot switch, exhaust blower, exhaust duct, drain valve, waterproof sockets*2, pipe strap*2			
Optional Accessory	Water and gas tap, Electric height adjustable base stand; Wind speed self adjustment function, real-time compensation for wind speed through sensor monitoring; APP networking remote control; USB data transfer function, supporting data monitoring and export; LED fluorescent lamp brightness adjustment function; Base bracket with manual or electric height adjustment function; Manual lifting front window function; Manual automatic lifting front window function.			
Package	Main Body	1250*980*1765mm	1640*980*1825mm	2000*1040*1895mm
	Exhaust Blower (W*D*H)	920*790*690mm	920*790*790mm	920*790*840mm
Gross Weight	265kg;50kg	325kg;56kg	420kg;60kg	

Version 2024 034/035