


EN Certified Biological Safety Cabinet



Color Touch Screen
Touch operation, real-time and clearer display of various values and appointment timing function.

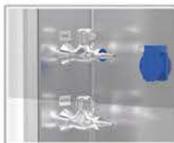

SS Water tap & Gas tap



Anti-paper dust structure




Wind Speed Sensor



Footmaster Caster
Universal caster with brake and leveling feet.




Stepper Motor

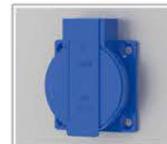

European Standard Socket



Work Zone
Work zone, made of 304 stainless steel, is surrounded by negative pressure.

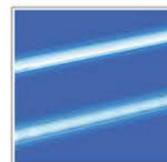

UV Lamp
Emission of 253.7 nanometers for most efficient decontamination.



ADVANTAGE

- 7-inch large touch color display: touch operation, real-time and clearer display of various values of the safety cabinet, dynamic air flow pattern diagram, fault alarm sound and light prompt.
- The power-on password can be set to prevent mis-operation by irrelevant personnel.
- Automatically adjust the wind speed.
- Power-off memory, restore the state before power-off, and have a visual reminder to alarm.
- Side windows on both sides for easy observation.
- The operation area is integrally formed, which is easy to disassemble and clean.
- A breeze speed sensor is equipped at the working area and the outer exhaust air outlet to detect the wind speed in real time.
- It can be used both manually and electrically, and the glass door can be pulled down to close the front window when the power is off.
- The cabinet is separated from the base. The height of the adjustable base can be customized.
- Interlocking function: UV lamp and front window; UV lamp and fan, LED lamp; fan and front window.
- It has an appointment timing function, which can automatically set the timing of turning on, turning off and the running time of the UV lamp, fan, and socket.
- The closing of the front window can double trigger the signal, so that the sterilization and disinfection function of the UV lamp can be turned on normally.
- The cabinet is designed with a 10° inclination angle, which conforms to the principle of ergonomics, has a larger viewing angle, is easy to operate and is more user-friendly.
- Anti-paper dust structure at the air inlet to prevent paper dust from entering the air cavity.
- Split countertop is easier to clean.
- One button front window lifting function.

TECHNICAL PARAMETERS

Model	BSC-3FA2-HA	BSC-4FA2-HA	
External Size(W*D*H)	1087*800*2260mm	1383*800*2260mm	
Internal Size(W*D*H)	980*600*650mm	1280*600*650mm	
Work Surface Height	750mm		
Max Opening	500mm		
Tested Opening	200mm		
Airflow Volume	Inflow	370m ³ /h(218cfm)	490m ³ /h(289cfm)
	Down flow	65%: 690 m ³ /h(406cfm)	65%: 910 m ³ /h(536cfm)
	Exhaust	35%: 370m ³ /h(218cfm)	35%: 490m ³ /h(289cfm)
Front Window	Manual and electric, glass thickness> 5mm, anti-ultraviolet		
ULPA Filter	2pcs, filtration efficiency for 0.12µm particles≥99.9995%		
Illumination	≥1000Lux		
Noise	EN12469≤65dB		
Display	Color Touch Screen		
European Standard Socket	230V, 3A		
Power Supply	AC 230V, 50/60Hz		
Rated Power	1100W	1200W	
Alarm	Alarm for abnormal wind speed, alarm for glass door not at safe height, alarm for abnormal power failure		
Motor	AC motor, 230V acceptable, adjustable speed, high efficiency, low power		
Material	Working Area: 304 stainless steel, side window glass; Frame and decorative plate: cold-rolled steel with anti-bacteria powder coating		
LED Lamp	12W*2	14W*2	
UV Lamp	20W	30W	
Standard Accessory	Cabinet*1, Base stand*1, LED lamp*2, UV lamp*2, S.S. water tap*1, S.S. gas tap*1, Drain ball valve*1, European Standard Socket*2, Remote control		
Optional Accessory	Arm rest, Formalin fumigation sterilizer, Infrared sterilizer, Electric height adjustable base stand, S.S. water tap*1, S.S. gas tap*1, Foot switch		
Package Size (W*D*H)	1240*1060*1880mm	1530*1060*1880mm	
Gross Weight	270kg	345kg	