



**NSF Certified  
Class II A2  
Biological Safety  
Cabinet**



Display for NA series




Display for GL series

**Color Touch Screen**  
Touch operation, real-time and clearer display of various values and appointment timing function.




**Exhaust Thimble canopy (Optional)**

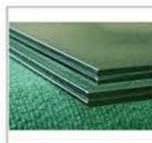

**Front Window**  
Two-layer laminated toughened glass >6mm, Anti UV



**Footmaster Caster**  
Universal caster with brake and leveling feet.

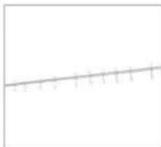

**Stainless Steel I.V. Bar with 10 hooks**



**SS Water Tap & Gas Tap**




**Waterproof Socket**



**Wind Speed Sensor**



**Work Zone**  
Unibody design

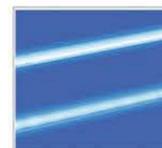

**UV Lamp**  
Emission of 253.7 nanometers for most efficient decontamination.



## ADVANTAGE

1. Motorized front window.
2. 7-inch large touch color display: touch operation, real-time and clearer display of various values of the safety cabinet, dynamic air flow pattern diagram, fault alarm sound and light prompt.
3. The power-on password can be set to prevent mis-operation by irrelevant personnel.
4. EC Fan, automatically adjust the wind speed.
5. With power-off memory function, the state before power-off will be restored after power-on and there will be visual prompts and alarms.
6. The integrated design of the operation area makes it easy to disassemble and clean.
7. The cabinet is separated from the base. The height of the adjustable base can be customized.
8. Interlocking function: UV lamp and front window; UV lamp and fan; fan and front window.
9. It has the function of appointment timing, which can automatically set the time of starting and shutting down the safety cabinet and the running time of ultraviolet lamps, fan and sockets.
10. The front window closes the trigger signal, so that the sterilization and disinfection function of the ultraviolet lamp is normally turned on.
11. The one-piece hand resting surface has a larger contact area with the arm and is more comfortable. air speed adjustable with filter block.

## TECHNICAL PARAMETERS

Model	BSC-4FA2	BSC-3FA2-GL BSC-3FA2-NA	BSC-4FA2-GL BSC-4FA2-NA	BSC-6FA2-GL BSC-6FA2-NA	
External Size(W*D*H)	1383*775*2220mm	1087*775*2220mm	1383*775*2220mm	1873*775*2240mm	
Internal Size(W*D*H)	1210*600*660 mm	910*600*660 mm	1210*600*660 mm	1700*600*660 mm	
Work Surface Height	762mm	762mm	762mm	782mm	
Max Opening	470mm				
Tested Opening	Safety height 200mm (8")		Safety height 254mm (10")		
Airflow Volume	Inflow	462m³/h (272cfm)	445m³/h(259cfm)	587m³/h(346cfm)	824m³/h(485cfm)
	Downflow	65%: 863m³/h (507cfm)	60%: 649m³/h(382cfm)	60%: 863m³/h(507cfm)	60%: 1212m³/h(713cfm)
	Exhaust	35%: 462m³/h (272cfm)	40%:445m³/h(259cfm)	40%: 587m³/h(346cfm)	40%: 824m³/h(485cfm)
Front Window	Motorized. Two-layer laminated toughened glass ≥6mm. Anti UV.				
ULPA Filter	Two, 99.9995% efficiency at 0.12µm, filter life indicator.				
Filter Guard Type	Aluminum alloy frame, glass fiber filter paper				
Noise	NSF 49s67dB				
Illumination	≥1000Lux				
Display	Color Touch Screen				
Waterproof Socket	Two, Total consumption: ≤500W				
Tap	Water tap*1, Gas tap *1				
Ground Resistance	≤ 0.10Ω				
Power Supply	AC230V, 50/60Hz; AC115V, 60Hz, Full load Amps: 9A, BTU/Hr: 1689				
Consumption	500W	400W	500W	700W	
Alarm	Abnormal airflow velocity; Filter replacement; Front window at unsafe height; High filter pressure alarm				
Motor	EC motor, 115V& 230V acceptable, Speed adjustable, high efficiency and low power consumption				
Material	Working area: 304stainless steel				
	Frame and Decorative Plate: cold-rolled steel with anti-bacteria power coating				
UV Lamp	30W	20W	30W	40W	
	Germicidal UV lamp, Emission of 253.7 nanometers for most efficient decontamination				
Illuminating Lamp	20W*2	10W*2	20W*2	23W*2	
	LED lamp				
Standard Accessory	LED lamp*2, UV lamp*2, Base stand, SS water tap*1, SS gas tap*1, Drain valve,				
	Waterproof sockets*2, Stainless Steel IV Bar with 10 hooks,Remote control				
Optional Accessory	Armrest, Airflow tester, Formalin fumigation sterilizer,Infrared sterilizer,				
	Electric height adjustable base stand				
Package Size (W*D*H)	1540*1080*1900mm	1240*1060*1900mm	1540*1080*1900mm	2020*1080*1900mm	
Gross Weight	330kg	280kg	330kg	410kg	

Version 2024 009/010